



# Surgical masks and caps



# Professional protection by evercare®

Surgical masks are intended to be worn by health care professionals during surgery. They should be worn at other times to catch the bacteria shed in liquid droplets and aerosols from the wearer's mouth and nose.

European Standard EN 14683 specifies the construction and performance requirements, and test methods for surgical masks intended to limit the transmission of infective agents from staff to patients and (in certain situations vice-versa) during surgical procedures in operating theatres and other medical settings with similar requirements.

## Key requirements of EN14683

### 1. Bacterial Filtration Efficiency

Determination of the filtration efficiency of the face mask. The higher the BFE %, the better the protection.

Classification: BFE =  $\geq 95$  % TYPE I  
BFE =  $\geq 98$  % TYPE II

### 2. Breathing Resistance

Determination of the airflow resistance of the surgical mask. The lower the differential pressure, the higher the breathing comfort.

Classification: TYPE I & II (non splash resistant) =  $< 29.4$  Pa/cm<sup>2</sup>  
TYPE IR & IIR (splash resistant) =  $< 49.0$  Pa/cm<sup>2</sup>

### 3. Splash Resistance

Determination of the resistance of the mask of potentially contaminated fluid splashes.

Classification: TYPE I & TYPE II not applicable  
TYPE IR & TYPE IIR  $\geq 120$  mmHg\*

\*120 mmHg is a minimum value. It corresponds to the average systolic arterial blood pressure, and intends to protect against ruptures in small arteries causing small sprays of blood.

evercare® surgical mask assortment includes a selection of Type II and Type IIR surgical masks.  
evercare® surgical masks complies to the EN14683:2005 for surgical masks.

evercare® surgical masks are soft, comfortable and provide excellent fit.



Surgical masks and caps  
Surgical masks



**Multiguard**

**REF 2850**

Type II  
Tie-band, blue  
Pack count: 50/300



**Multifix**

**REF 2856**

Type II  
Ear-loop, green  
Pack count: 50/300



**Multiguard**

**REF 2851**

Type II  
Tie-band, green  
Pack count: 50/300



**Multifix Guardline**

**REF 2858**

Type II  
Ear-loop, green  
Pack count: 50/300



**Multiguard Antifog**

**REF 2865**

Type II  
Tie-band, antifog, green  
Pack count: 50/300



**Multifix Splash**

**REF 2875**

Type IIR  
Ear-loop, splash resistant,  
green  
Pack count: 50/300



**Multifix Starline**

**REF 2855**

Type II  
Ear-loop, blue  
Pack count: 50/300



**Duckbill**

**REF 2852**

Type II and Type IIR  
Tie-band, splash-resistant,  
Duckbill, blue/white  
Pack count: 50/300



# Surgical masks & hoods

Surgical masks and caps



## Surgical Mask

**REF 2853**

Type IIR  
Tie-band, extra wide,  
splash-resistant, anti-fog,  
soft, white  
Pack count: 50/300



## Surgical hood Gallant, shoulder protection

Patterned viscose  
nonwoven

**REF 2779**

Green  
Pack count: 6 x 50



## Surgical Mask Face shield

**REF 2854**

Type IIR  
Tie-band, visor, anti-fog,  
anti-reflection, soft,  
splash resistant, blue  
Pack count: 25/100



## Surgical hood Lancelot, neck protection

Patterned viscose  
nonwoven

**REF 2785**

Green  
Pack count: 6 x 50



## Surgical Mask

**REF 2876**

Type IIR  
Ear-loop, splash resistant,  
soft, anti-fog, white  
Pack count: 50/300



## Surgical hood James, tie-on

Patterned viscose  
nonwoven

**REF 2773**

Green  
Pack count: 6 x 100



## Surgical Mask, Soft, type II

**REF 2860**

Type II  
Tie-band, extra wide,  
soft, anti-fog, white  
Pack count: 50/300



## Surgical hood James, tie-on

Patterned viscose  
nonwoven

**REF 2774**

Blue  
Pack count: 6 x 100

# Surgical hoods and caps



Surgical hood  
Valiant, shoulder  
protection

Patterned viscose  
nonwoven

**REF 2780**

Size: Universal  
Blue  
Pack count: 6 x 50

**REF 2796**

Size: S  
Blue  
Pack count: 6 x 50



**REF 2781**

Size: XL  
Blue  
Pack count: 6 x 50

**REF 2797**

Size: Universal  
Violet  
Pack count: 6 x 50



Surgical cap  
Tricolor

Patterned viscose  
nonwoven, soft  
spunlace sweatband

**REF 2770**

Green/blue/white  
Pack count: 6 x 50



Surgical cap  
Kolumbus

Patterned viscose  
nonwoven

**REF 2776**

Green  
Pack count: 6 x 100



Surgical cap  
Lady

Patterned viscose  
nonwoven

**REF 2794-01**

Green  
Pack count: 6 x 100



Surgical cap  
Kolumbus

Patterned viscose  
nonwoven

**REF 2777**

Blue  
Pack count: 6 x 100



Surgical cap  
Julia,  
for long hair

Patterned viscose  
nonwoven

**REF 2784**

Violet  
Pack count: 6 x 50



Surgical cap  
Kolumbus

Patterned viscose  
nonwoven

**REF 2778**

Violet  
Pack count: 6 x 100

Surgical masks and caps

# Surgical caps



## Surgical cap Romeo

Patterned viscose nonwoven

**REF 2788**

Green  
Pack count: 6 x 100



## Round cap Rondo Standard

Polypropylene nonwoven

Pack count: 6 x 150

**REF 2740-01**

M Green

**REF 2741-01**

M White

**REF 2742-01**

M Blue

**REF 2743-01**

M Pink

**REF 2744-01**

M Multicolor  
(Green/Blue/Pink)

**REF 2750-01**

L Green

**REF 2751-01**

L White



## Surgical cap Tristan, tie-on

Patterned viscose nonwoven

**REF 2791**

Green  
Pack count: 6 x 100



**REF 2752-01**

L Blue

**REF 2753-01**

L Pink

**REF 2754-01**

L Multicolor  
(Green/Blue/Pink)

**REF 2760-01**

XL Green



## Surgical cap Tina

Patterned viscose nonwoven

**REF 2782**

Green  
Pack count: 6 x 100



**REF 2761-01**

XL White

**REF 2762-01**

XL Blue

**REF 2763-01**

XL Pink

**REF 2764-01**

XL Multicolor  
(Green/Blue/Pink)



## Surgical cap Tina

Patterned viscose nonwoven

**REF 2783**

Blue  
Pack count: 6 x 100



# Surgical caps & Dispenser holders



## Rondo Special

Polypropylene nonwoven

**REF 2765**

L Blue  
Pack count: 6 x 150



## Rondo Extra

Polypropylene nonwoven

**REF 2766**

L Blue  
Pack count: 6 x 100


## Dispenser holders for Surgical Masks and Headwear

### Features

- Adjustable to fit dispensers in the evercare® mask and headwear assortment
- Available in 2 sizes; “small dispensers and large dispensers”
- Made of stainless steel that enable you to clean the holder easily
- Convenient in the right place in tidy appearance
- Simplicity in refilling; just change the box
- Easy to access the mask or headwear from the dispenser

### Use

- Outside the operating room, on the wall, next to the disinfection/hand scrub area
- At the entrance to the Intensive Care Unit
- In emergency/acute care centers examination rooms
- In the microbiological laboratories the sterilization departments

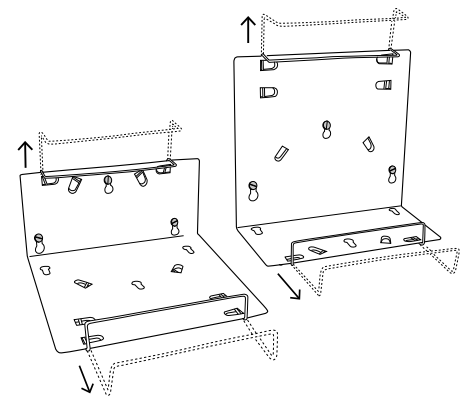

Surgical masks and caps



### Small

**REF 2800**

Pack count: 1 x 28

Suits masks:

2850, 2851, 2853, 2855, 2856, 2858, 2857, 2860, 2865, 2875, 2876

The grips are flexible and they hold the dispenser box firmly with tensioned wires.



### Large

**REF 2700**

Pack count: 1 x 8

Suits masks: 2852 and 2854

Suits caps:

Valiant universal and Valiant XL, Lancelot, Gallant, James, Kolumbus, Romeo, Julia, Tricolor, Tristan, Tina



## Label information in accordance with EN ISO 15223-1 and EN 1041

